



# AN 412

## Accuracy in SiON Thickness and N Dose Measurements by XPS

May 7, 2007 (Version 3.0)

### Summary

In summary, a series of comparisons have been performed between the SiON High Precision XPS Measurement and other analytical methods for both SiON thickness and N dose. The agreement between methods among all of the comparisons is good. Although none of the methods described can claim to give the “absolute” correct values, these results offer confidence that the SiON High Precision XPS analysis provides both measurements using a single technique while obtaining results comparable to those from all other techniques applicable for measuring these characteristics of SiON gate dielectric films.

### Discussion

SiON gate dielectric films for 130 nm, 90 nm and 65 nm node IC devices must be tightly controlled in both thickness and nitrogen dose, because these variables can not only be independent from one another, but each alone can affect the device performance. The SiON High Precision XPS Measurement has been developed to address these requirements, i.e., independent, high precision measurements of thickness and nitrogen dose in SiON gate dielectric films.

However, the accuracy of an analytical measurement is sometimes a point of discussion, especially when NIST traceable reference materials are unavailable. In this instance a comparison to other analytical methods can be helpful to build confidence in the measurement. The SiON High Precision XPS Measurement for SiON thickness and N dose is compared below to ellipsometry, TEM, NRA (Nuclear Reaction Analysis) and SIMS.

For thickness comparison to ellipsometry, thin SiO<sub>2</sub> films were used since ellipsometry measurement of SiON films is unreliable. Table 1 (upper half) shows the comparison for three samples with thicknesses ranging from 5 to 24 Å. The agreement is good between the two techniques.

For thickness comparison to TEM, SiON films were used. The TEM cross sections for four SiON films are shown in Figures 1 to 4. The comparison of the thickness values between TEM and XPS is shown in Table 1 (lower half) for thicknesses ranging from 17 to 56 Angstroms. The thickness values are consistent between the two techniques. Note that the TEM thickness values depend upon where the “interface” and “top surface” of the SiON film are defined (see Figures 2 and 4). For a 15Å SiON film the distance between these two planes can vary by 10-30% across the TEM section image. In contrast the SiON High Precision XPS Measurement of SiON thickness has a precision <1% RSD.

**Table 1. Thickness Comparison: High Precision XPS, Ellipsometry and TE**

	SiO <sub>2</sub> High Precision XPS and Ellipsometry	
	XPS (wafer center)	Ellipsometry (nominal wafer avg.)
SiO <sub>2</sub>		
Wafer A	5.55 Å	6.0 Å
Wafer B	15.66 Å	16.1 Å
Wafer C	24.30 Å	24.5 Å
SiON High Precision XPS and TEM		
SiON	XPS	TEM
Wafer D	17.40 Å	18 Å
Wafer E	24.74 Å	27 Å
Wafer F @ 0.10mm	46.83 Å	43 Å
Wafer F @ 0.90mm	56.35 Å	54 Å

For a N dose comparison to NRA, four SiON samples with N doses ranging from  $6E14$  atoms/cm<sup>2</sup> to  $2.5E15$  atoms/cm<sup>2</sup> were used. The results are shown in Table 2. Again, good agreement exists between the two techniques.

For a N dose comparison to SIMS, a low dose nitrogen implant into SiO<sub>2</sub> was used. The result is shown in Table 3. Good consistency is found between SIMS and the XPS measurement for N dose.

In summary, a series of comparisons have been performed between the SiON High Precision XPS Measurement and other analytical methods for both SiON thickness and N dose. The agreement between methods among all of the comparisons is good. Although none of the methods described can claim to give the “absolute” correct values, these results offer confidence that the SiON High Precision XPS analysis provides both measurements using a single technique while obtaining results comparable to those from all other techniques applicable for measuring these characteristics of SiON gate dielectric films.


**Table 2. N Dose Comparison: SiON High Precision XPS and NRA**

	Thickness (Å)	XPS (atoms/cm <sup>2</sup> )	NRA (atoms/cm <sup>2</sup> )	NRA Precision
Wafer D	17.33	6.37E+14	6.90E+14	(+/-0.69E+14)
Wafer G	26.39	2.99E+15	2.42E+15	(+/-0.24E+15)
Wafer E	24.74	5.86E+15	4.77E+15	(+/-0.48E+15)
Wafer H	18.73	2.53E+15	2.04E+15	(+/-0.20E+15)


**Table 3. N Dose Comparison: SiON**

	Thickness (Å)	XPS (atoms/cm <sup>2</sup> )	SIMS (atoms/cm <sup>2</sup> )
Wafer I	43.38	3.57E+14	3.70E+14


#### TEM Cross Section

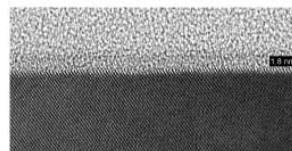


**Figure 1. Wafer D**

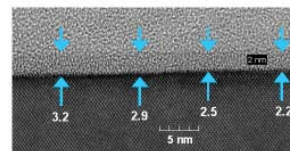

**Figure 2. Wafer E**

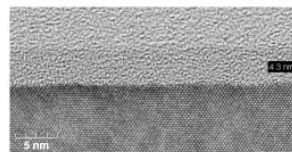

**Figure 3. Wafer F @ 0,10mm**

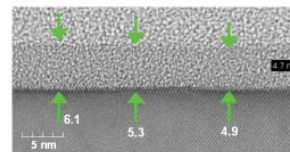

**Figure 4. Wafer F @ 0,90mm**

#### United States Locations

Tempe, Arizona  
+1 480 239 0602 info.az@eaglabs.com  
+1 602 470 2655 fax

Sunnyvale, California  
810 Kifer Road  
+1 408 530 3500 info.ca@eaglabs.com  
+1 408 530 3501 fax

1135 E Arques Avenue  
+1 408 738 3033  
+1 408 738 3035 fax

785 Lucerne Drive  
+1 408 737 3892  
+1 408 737 3916 fax

Peabody, Massachusetts  
+1 978 278 9500 info.ma@eaglabs.com  
+1 978 278 9501 fax

Chanhassen, Minnesota  
+1 952 828 6411 info.mn@eaglabs.com  
+1 952 828 6449 fax

East Windsor, New Jersey  
+1 609 371 4800 info.nj@eaglabs.com  
+1 609 371 5666 fax

Syracuse, New York  
+1 315 431 9900 info.ny@eaglabs.com  
+1 315 431 9800 fax

Raleigh, North Carolina  
+1 919 829 7041 info.nc@eaglabs.com  
+1 919 829 5518 fax

Round Rock, Texas  
+1 512 671 9500 info.tx@eaglabs.com  
+1 512 671 9501 fax

#### International Locations

Shanghai, China  
+ 86 21 6879 6088 info.cn@eaglabs.com  
+ 86 21 6879 9086 fax

Tournefeuille, France  
+ 33 5 61 73 15 29 info.fr@eaglabs.com  
+ 33 5 61 73 15 67 fax

Frankfurt, Germany  
+ 49 (0) 693053213 info.de@eaglabs.com  
+ 49 (0) 69307941 fax

Tokyo, Japan  
+ 81 3 5396 0531 info.jp@eaglabs.com  
+ 81 3 5396 1930 fax

HsinChu, Taiwan  
+ 886 3 5632303 info.tw@eaglabs.com  
+ 886 3 5632306 fax

Uxbridge, United Kingdom  
+ 44 (0) 1895 811194 info.uk@eaglabs.com  
+ 44 (0) 1895 810350 fax